


RenuLife Violet Ray

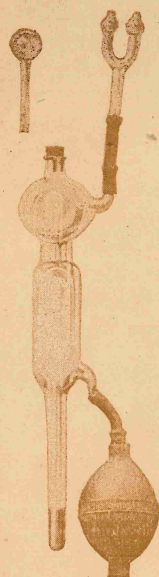

No. 25 Ozone Generator

Ozone is simply a concentrated form of Oxygen. When we breathe, we fill the lungs full of air or oxygen. This oxygen passes through the walls of the lungs, at the same time passing through and purifying the blood, which is flowing through the walls of the lungs and which comes there from every part of the body to be purified. It is the lack of pure air to purify our blood that causes us to be overcome with disease. Being forced to spend many hours daily within the confines of an office or factory, is responsible for most of our diseases, and particularly diseases of the respiratory organs. The RenuLife Ozone Generators deliver Ozone in large quantities for inhalation purposes. The Ozone penetrates every cell in the lungs, purifying the blood, destroying disease germs and soothing the inflamed tissues. For ASTHMA, CATARRH, HAY FEVER, ROSE FEVER, COUGHS, COLDS, BRONCHITIS, etc., there is no better treatment known. Our Ozone Inhaler attachment No. 25 is designed so as to fit the handle of any of our Generators. The High Frequency current conducted through the attachment creates Ozone in quantities, confines it, and filters it through an oil so that one may inhale it with maximum benefits.

Price.....\$9.50

ELECTRODES FOR INTERNAL TREATMENTS

For internal treatment, electrodes are furnished to fit any orifice of the body. We furnish them in plain or the insulated type. The insulated electrode has a small vacuum tube running through the center of the electrode, and all current is brought to the tip of the electrode, directly where the current is desired. Plain electrodes deliver the current evenly over all parts of the electrode touched by the body during treatment. We recommend that all internal electrodes should be insulated. Best results will be obtained from them. See Nos. 5, 6, 7, 8, 9, 10, 12, 21, 22, 23, 24 and 27.

Dental Electrodes

Dental electrodes are supplied in five shapes, and are particularly designed for the treatment of the gums and teeth. Electrode No. 40 is for cavity treatment, and for diagnosing live or dead nerves—No. 41 is for the treatment of the frontal gums—No. 42 is a cataphoric electrode for carrying medicaments in the treatment of abscesses of the back gums—No. 43 is used for the same purpose as No. 42, except it is for treatment of the frontal gums—No. 44 is a special PYORRHEA electrode, and is used by hundreds of dentists throughout the country with splendid results in the treatment of PYORRHEA. It restores circulation to the gums, stopping the pus, and restoring circulation of the gums to normal. By doing this, the gums take on a healthy appearance, and fill in where they have started to recede.

Price of Each Dental Electrode.....\$2.50

High Frequency Liniment

This preparation is manufactured by the RenuLife Electric Co. especially for use in treatment of RHEUMATISM, NEURALGIA, NEURITIS, LUMBAGO, etc. We particularly recommend this treatment for use in connection with the RenuLife Violet Ray. A small quantity of it rubbed on the body, followed by a Violet Ray treatment, will be driven directly into the tissues, with beneficial results.

Price, per Bottle.....\$1.00

Doctor Eberhart's Manual—A Working Manual of High Frequency Currents

By Noble M. Eberhart, M. D., Ph. D., D.C.L., Professor and Head of Department of Physiologic Therapeutics, Med. Dept., Loyola University, Chicago; Professor of High Frequency, Illinois School of Electro-Therapeutics.

Just what the name implies—a working handbook, telling how to use these currents. No lengthy discussions, no perplexing technicalities (no red tape). 324 pages, over 60 pages of illustrations. Handsomely bound.

Price.....\$4.00

RenuLife Violet Ray

GENERAL INSTRUCTIONS FOR THE USE AND CARE OF RENU LIFE ELECTRODES

In treating the surface of the body, always keep the electrode moving slowly. Never insert or remove an electrode from an orifice of the body when the current is turned on. Electrodes that are to be used internally should be lubricated with cold cream or vaseline before being inserted. High Frequency electricity is a germicide, yet sanitary rules must be followed in the care of electrodes. To sterilize them, immerse in a weak solution of carbolic acid or in Lysol. Then wipe thoroughly with a clean cloth or absorbent cotton. IMPORTANT—To get the maximum strength from an electrode, always hold the end of handle nearest cord.

Always insert electrodes into the handle with a straight push. Don't TWIST. Remove with a straight pull. DON'T TWIST. Never use the generator when in a bath tub or when grounded by means of a radiator, electric fixtures, water pipes, etc.

Never use the generator on the scalp when the hair has been treated with a tonic containing more than 40% alcohol.

Should the electrode at any time not give off a spark when the generator is operating, the fault is with the electrode. Air has gotten into it and destroyed the vacuum. It has very likely in some way been cracked, although the crack may not be visible. In any event, if the generator or electrodes are not working properly, we would suggest that you write us immediately, telling us just exactly what the conditions are. This may avoid the necessity of returning generator to us.

Talcum Powder should be dusted over the surface of the body being treated to allow free movement of glass electrodes.

INSTRUCTIONS FOR OPERATING OUR OZONE GENERATOR No. 25

Ozone Generators are packed with nose piece attached: included separately is a small glass piece with flat disc end for inhalations through mouth. The nose piece may be removed from rubber coupling and mouthpiece inserted. Place rubber bulb on glass tip at lower part of barrel of generator. Place one teaspoonful of Pine Needle Oil in upper chamber by removing cork. This oil will last for a considerable length of time.

Insert stem of Ozone Generator in handle of Violet Ray Generator, same as any other electrode. Hold barrel chamber of Ozone Generator firmly in left hand. This causes sparks to jump from inner vacuum to outer wall, creating the Ozone. Turn on current and agitate bulb with right hand. This forces the Ozone into the upper chamber, where it is filtered through the oil. A pleasing, purifying, powerful agency results.

INSTRUCTIONS FOR USING TONIC SATURATOR ELECTRODE

Turn on current by means of adjusting knob. Grasp electrode firmly in the hand and insert in handle of generator. Pull out electrode before turning off current. Do not use this electrode when in contact with metal or when grounded by means of water or electric fixtures. Do not attempt to adjust controlling knob with one hand before pulling the electrode from the handle.

Electrode Applicators

For Use in Connection With

RenuLife Violet Ray Generators

1304 Harper Ave
Madison 2940

Manufactured by

Beasley-Eastman Laboratories
INCORPORATED


4841-53 Rivard Street
DETROIT, MICHIGAN

Originators and Manufacturers of
RenuLife Violet Ray Generators

RENULIFE VIOLET RAY SALES CO.,
512 EMPIRE BUILDING,

OLIFFORD AT WASHINGTON BLVD.,

DETROIT, MICH.

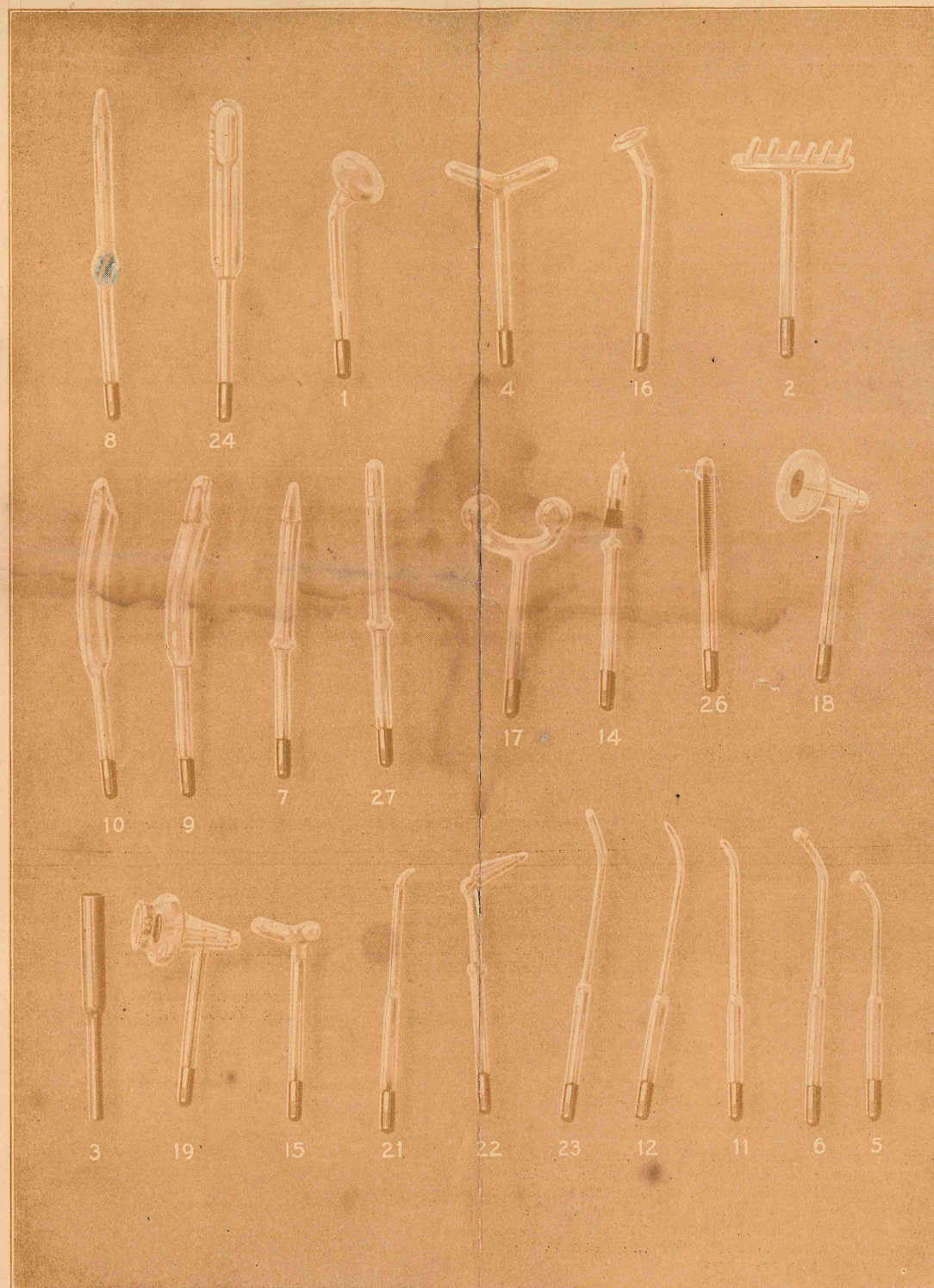

- No. 1 Surface or General Body Electrode**
Suited for all exterior surface treatments.....\$1.00
- No. 2 Rake Comb Electrode**
Especially for scalp treatment, Falling Hair, Dan-
druff..... 2.25
- No. 3 Metal Electrode**
Held tightly in hand for systematic treatment..... 1.00
- No. 4 External Throat Electrode**
Specially designed to fit closely over neck, for Goi-
tre, Tonsilitis..... 1.75
- No. 5 Internal Throat Insulated Electrode**
Particularly for professional use..... 2.00
- No. 6 Internal Throat Plain Electrode**
Particularly for professional use..... 1.25
- No. 7 Rectal Insulated Electrode**
Designed for treatment of prostatic trouble..... 2.00
- No. 8 Rectal Plain Electrode**
Also for rectal and prostatic treatment, professional
use..... 1.25
- No. 9 Vaginal Insulated Electrode**..... 2.00
- No. 10 Vaginal Plain Electrode**..... 1.25
- No. 11—Discontinued**
- No. 12 Urethral Plain Electrode**
For professional use only..... 2.00
- No. 14 Fulguration Electrode**
For cauterizing and removing blemishes, scars,
warts, etc..... 1.50
- No. 15 Spinal Electrode**
Best for spinal treatments..... 1.75
- No. 16 Single Eye Electrode**
Made to fit over eye, with lid closed..... 1.25
- No. 17 Double Eye Electrode**
For treating both eyes at the same time..... 2.00
- No. 18 Condenser Electrode**
Made with metal conductor within vacuum to con-
centrate current, producing strong, even flow. Best
for deep seated cases in place of No. 1..... 3.50
- No. 19 Cataphoric Electrode**
For professional use..... 3.50
- No. 21 Ear and Nasal Insulated Electrode**
For internal ear treatment for deafness; also for in-
ternal nose treatment..... 2.00
- No. 22 Ear Special Insulated Electrode**
Best suited for internal ear treatments..... 2.25
- No. 23 Ear and Nasal Plain Electrode**
For internal treatment of nose and ear, for profes-
sional use..... 1.25
- No. 24 Special Vaginal Electrode**
Made with perforations, generating purifying ozone
internally. Particularly for inflammation..... 2.75
- No. 25 Ozone Generator**
See next page for description..... 9.50
- No. 26 RenuLife Tonic Saturator Electrode**
Specially designed for saturation treatment..... 2.50
- No. 27 Special Prostatic Electrode**..... 2.50