

No. 25 Ozone Generator

Ozone is simply a concentrated form of Oxygen. The RenuLife Ozone Generators deliver Ozone in large quantities for inhalation purposes.

Our Ozone Inhaler attachment No. 25 is designed so as to fit the handle of any of our Generators. The High Frequency current conducted through the attachment creates Ozone in quantities, confines it, and filters it through an oil so that one may inhale it.

Price\$12.50

INSTRUCTIONS FOR OPERATING OUR OZONE GENERATOR No. 25

Ozone Generators are packed with nose piece attached: included separately is a small glass piece with flat disc end for inhalations through mouth. The nose piece may be removed from rubber coupling and mouthpiece inserted. Place rubber bulb on glass tip at lower part of barrel of generator. Place one teaspoonful of Pine Needle Oil in upper chamber by removing cork. This oil will last for a considerable length of time.

Insert stem of Ozone Generator in handle of Violet Ray Generator, same as any other electrode. Hold barrel chamber of Ozone Generator firmly in left hand. This causes sparks to jump from inner vacuum to outer wall, creating the Ozone. Turn on current and agitate bulb with right hand. This forces the Ozone into the upper chamber, where it is filtered through the oil.

ELECTRODES FOR INTERNAL TREATMENTS

For internal treatment, electrodes are furnished to fit any orifice of the body. We furnish them in plain or the insulated type. The insulated electrode has a small vacuum tube running through the center of the electrode, and all current is brought to the tip of the electrode directly where the current is desired. Plain electrodes deliver the current evenly over all parts of the electrode touched by the body during treatment. We recommend that all internal electrodes should be insulated. Best results will be obtained from them. See Nos. 5, 6, 7, 8, 9, 10, 12, 21, 22, 23, 24 and 27.

High Frequency Liniment

This preparation is manufactured by the RenuLife Electric Co. especially for use in treatment of RHEUMATISM, NEURALGIA, NEURITIS, LUMBAGO, etc. We particularly recommend this treatment for use in connection with the RenuLife Violet Ray. A small quantity of it rubbed on the body, followed by a Violet Ray treatment, will be driven directly into the tissues, with beneficial results.

Price, per Bottle\$1.00

GENERAL INSTRUCTIONS FOR THE USE AND CARE OF RENULIFE ELECTRODES

In treating the surface of the body, always keep the electrode moving slowly. Never insert or remove an electrode from an orifice of the body when the current is turned on. Electrodes that are to be used internally should be lubricated with cold cream or vaseline before being inserted. High Frequency electricity is a germicide, yet sanitary rules must be followed in the care of electrodes. To sterilize them, immerse in a weak solution of carbolic acid or in Lysol. Then wipe thoroughly with a clean cloth or absorbent cotton. **IMPORTANT**—To get the maximum strength from an electrode, always hold the end of the handle nearest cord.

Always insert electrodes into the handle with a straight push. Don't TWIST. Remove with a straight pull. DON'T TWIST. Never use the generator when in a bath tub or when grounded by means of a radiator, electric fixtures, water pipes, etc.

Never use the generator on the scalp when the hair has been treated with a tonic containing more than 10% alcohol.

Should the electrode at any time not give off a spark when the generator is operating, the fault is with the electrode. Air has gotten into it and destroyed the vacuum. It has very likely in some way been cracked, although the crack may not be visible. In any event, if the generator or electrodes are not working properly, we would suggest that you write us immediately, telling us just exactly what the conditions are. This may avoid the necessity of returning generator to us.

Talcum Powder should be dusted over the surface of the body being treated to allow free movement of glass electrodes.

Electrode Applicators

For Use in Connection With

RenuLife Violet Ray Generators

Manufactured by

RenuLife Electric Co.

Harper at Rivard


DETROIT 11, MICHIGAN

*Originators and Manufacturers of
RenuLife Violet Ray Generators*

*Please include in your remittance
postage and sales tax (if applicable)
when ordering electrodes.*

*SEND EXPRESS OR P. O. MONEY ORDER
CURRENCY SENT AT YOUR RISK.*

Rec. April 19th 48

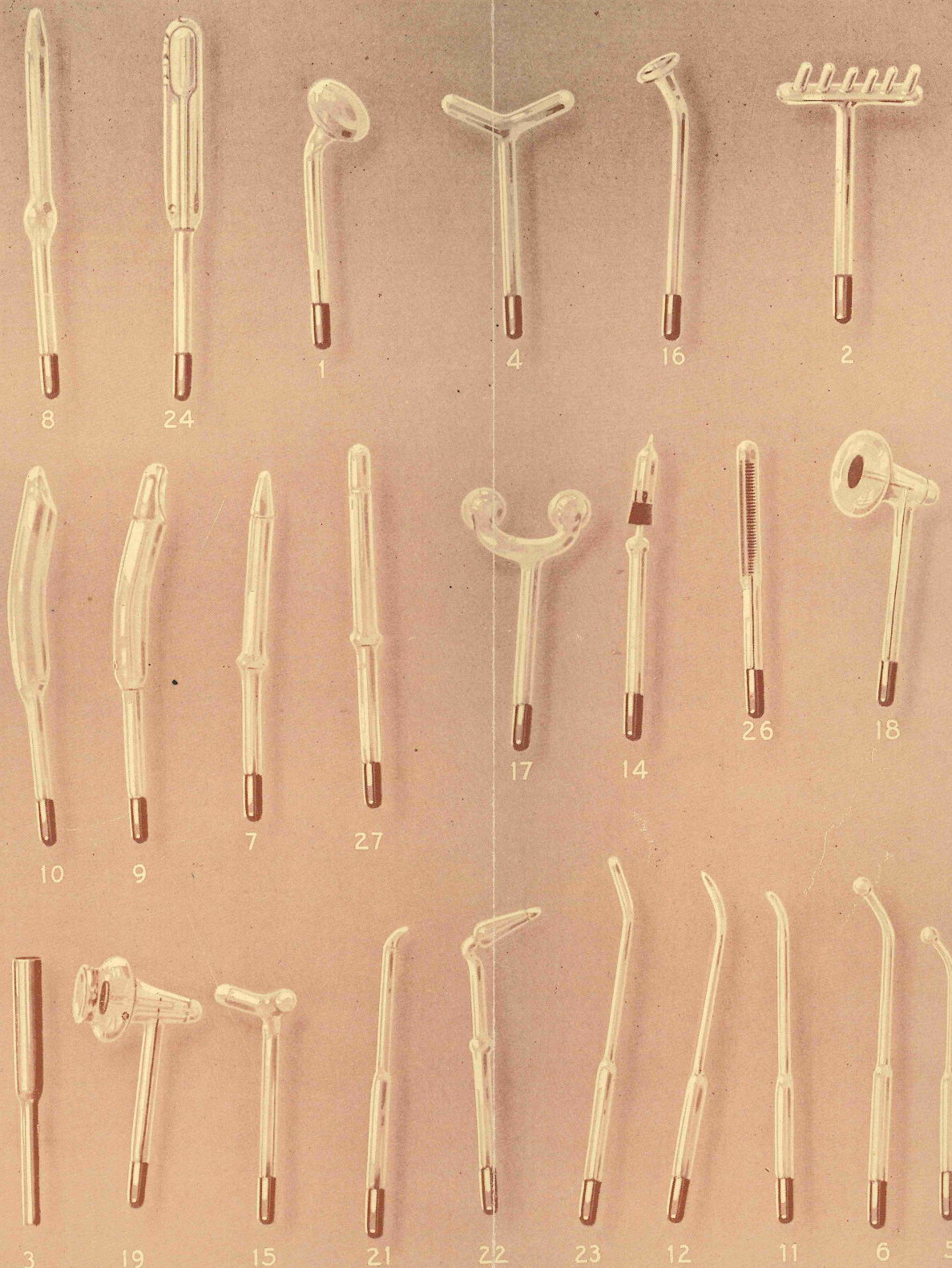

- No. 1 Surface or General Body Electrode**
Suited for all exterior surface treatments.....\$1.25
- No. 2 Rake Comb Electrode**
Especially for scalp treatment, Falling Hair, Dandruff..... 2.50
- No. 3 Metal Electrode**
Held tightly in hand for systemic treatment..... 1.50
- No. 4 External Throat Electrode**
Specially designed to fit closely over neck, for Goitre, Tonsillitis..... 1.75
- No. 5 Internal Throat Insulated Electrode**
Particularly for professional use..... 2.50
- No. 6 Internal Throat Plain Electrode**
Particularly for professional use..... 1.25
- No. 7 Rectal Insulated, Electrode**
Designed for treatment of prostratic trouble.... 3.50
- No. 8 Rectal Plain Electrode**
Also for rectal and prostatic treatment, professional use..... 1.75
- No. 9 Vaginal Insulated Electrode** 3.50
- No. 10 Vaginal Plain Electrode** 1.75
- No. 11 Discontinued**
- No. 12 Urethral Plain Electrode**
For professional use only..... 2.00
- No. 14 Fulguation Electrode**
For cauterizing and removing blemishes, scars, warts, etc..... 1.75
- No. 15 Spinal Electrode**
Best for spinal treatments..... 2.00
- No. 16 Single Eye Electrode**
Made to fit over eye, with lid closed..... 1.50
- No. 17 Double Eye Electrode**
For treating both eyes at the same time..... 2.75
- No. 18 Condenser Electrode**
Made with metal conductor within vacuum to concentrate current, producing strong, even flow. Best for deep seated cases in place of No. 1 4.50
- No. 19 Cataphoric Electrode**
For professional use..... 5.50
- No. 21 Ear and Nasal Insulated Electrode**
For internal ear treatment for deafness; also for internal nose treatment..... 2.75
- No. 22 Ear Special Insulated Electrode**
Best suited for internal ear treatments..... 3.50
- No. 23 Ear and Nasal Plain Electrode**
For internal treatment of nose and ear, for professional use..... 2.25
- No. 24 Special Vaginal Electrode**
Made with perforations, generating purifying ozone internally. Particularly for inflammation.. 3.50
- No. 25 Ozone Generator**
See next page for description.....12.50
- No. 26 Renulife Tonic Saturator Electrode**
Specially designed for saturation treatment.... 2.50
- No. 27 Special Prostrate Electrode** 3.50